
*Escali*

**PRODUCT  
MANUAL**

# **Glass Waves Body Scale**

Model # B180SG

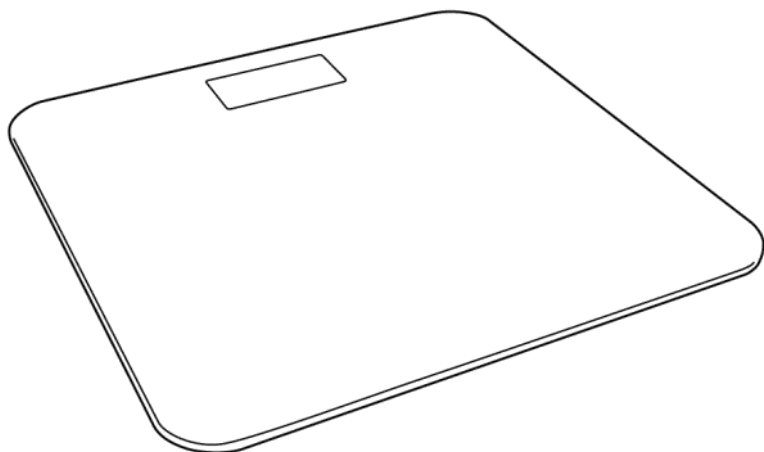**[www.escali.com](http://www.escali.com)**

**Thank you for choosing the Escali Glass Waves Body Scale.** This precision state-of-the-art measuring instrument is most accurate when weighing on a flat hard surface. This scale will not function correctly on carpet. As the electronic sensors are sensitive, avoid dropping, or jarring the scale, and store where it will be protected from impact.

## **OPERATING INSTRUCTIONS**

1. Remove the insulating plastic strip from the battery, and ensure the **+/-** terminals are correctly placed as shown inside the battery compartment.
2. Select your desired mode of display by pressing the **'unit'** button on the base of the scale.
3. Place the scale on a solid, and flat surface.
4. To activate, simply stand directly onto the scale platform.
5. The scale will turn off automatically when finished.

## **PROPER USE & CARE**

- Use and stow the scale at room temperature.
- Use on a flat, level, and stable surface.
- This scale is an instrument of precision. Do not submit the scale to shaking, excessive vibration, or other rough treatment.
- The display may be affected by electromagnetic disturbances such as radios, microwaves, or cell phones. If such disturbances occur, remove the source of the disturbance, and restart the scale.
- Disassembling, or tampering with the scale voids the warranty.
- Clean the scale with a damp cloth.
- Do not submerge in liquid. The scale is not waterproof.

- For best results, use lithium cell batteries.
- Install only new batteries of the same type in your product.
- Failure to insert batteries in the correct polarity, as indicated in the battery compartment, may shorten the life of the batteries or cause batteries to leak.
- Do not mix old and new batteries.
- Do not mix Alkaline, Standard (Carbon-Zinc) or Rechargeable (Nickel Cadmium) or (Nickel Metal Hydride) batteries.
- Do not dispose of batteries in fire.
- Batteries should be recycled or disposed of as per state and local guidelines.

## **WARNING MESSAGES**

Lo Low battery indicator

Err Overload indicator

## **PRODUCT SPECIFICATIONS**

Model: B180SG

Capacity: 400 lb / 180 kg

Graduation: 0.2 lb / 0.1 kg / 0.2 st

Product Dimensions: 11.9 x 11.9 x 0.7 inches

Power Source: 1 - CR2032 lithium battery (included)

## **ESCALI CORPORATION**

3203 Corporate Center Drive, Suite 150  
Burnsville, MN 55306  
USA

### **PRODUCT SUPPORT**

24/7 Online support: <https://escali.zendesk.com>

Email support: [feedback@escali.com](mailto:feedback@escali.com)

Phone support: 1-800-467-6408

### **LIFETIME LTD WARRANTY**

This scale is guaranteed to be free of defects in materials and workmanship. For more details visit: [www.escali.com/warranties](http://www.escali.com/warranties)

Escali's liability is limited to two (2) times the cost of the product.

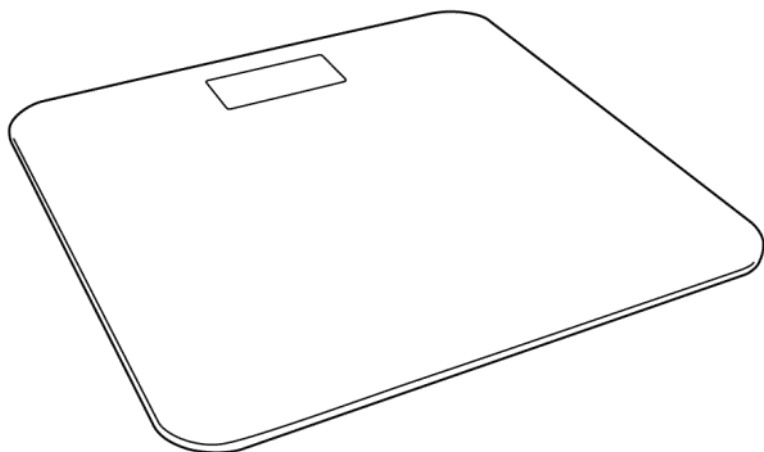
© Escali Corp., Minneapolis, MN, USA. All rights reserved. Escali, the Escali logo and other Escali marks are owned by Escali Corp. and may be registered.

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